

VERSAs FOOT₂



2021 vs2 Product Line Catalog
BIODAPTINC.COM | 320.252.0894



For competitive athletes and those weekend warriors who want to get the most out of their favorite sport or activities.

**REFUSE
TO COMPROMISE**



VERSA FOOT²



THE NEXT GENERATION VERSA FOOT

A versatile and highly adjustable prosthetic foot used for sport, recreational, or every day activities for both AK and BK amputees.

Purpose built to perform in the most extreme activities and environments, the design of the VF2 allows it to efficiently absorb shock and impact through the toe and heel.

Featuring customized sole plates for different activities and a wide range of adjustability and tuning options, the Versa Foot 2's unique ability to adapt makes it one of the most versatile feet available.

KEY FEATURES

Check out the significant performance improvements to the original Versa Foot, the choice of many competitive athletes including medal winning Paralympians.



4-BAR ANKLE LINKAGE LOWERS SHOCK LOCATION

Improves shock efficiency and lowers the overall profile allowing for a better fit inside a shoe or boot.



NEW HIGH PERFORMANCE FOX DPS FLOAT SHOCK

New tune and increased size.



NEW HEEL CUSHION

Allows vertical compression at heel strike.



HEEL SLIDER/GUARD

- Doubles as protective shield for shock adjuster.
- Allows the Versa Foot 2 to easily slide into a shoe or boot.




INTERCHANGEABLE CARBON FIBER SOLE PLATE

- 4 Options: Foam Sole, Shoe Sole, Bike Sole, MX Sole
- New, Vibram sole designs with increased durability to hold up better in the most extreme environments



THE NEXT GENERATION VERSA FOOT WITH INCREASED VERSATILITY.

HIGH PERFORMANCE FEATURES

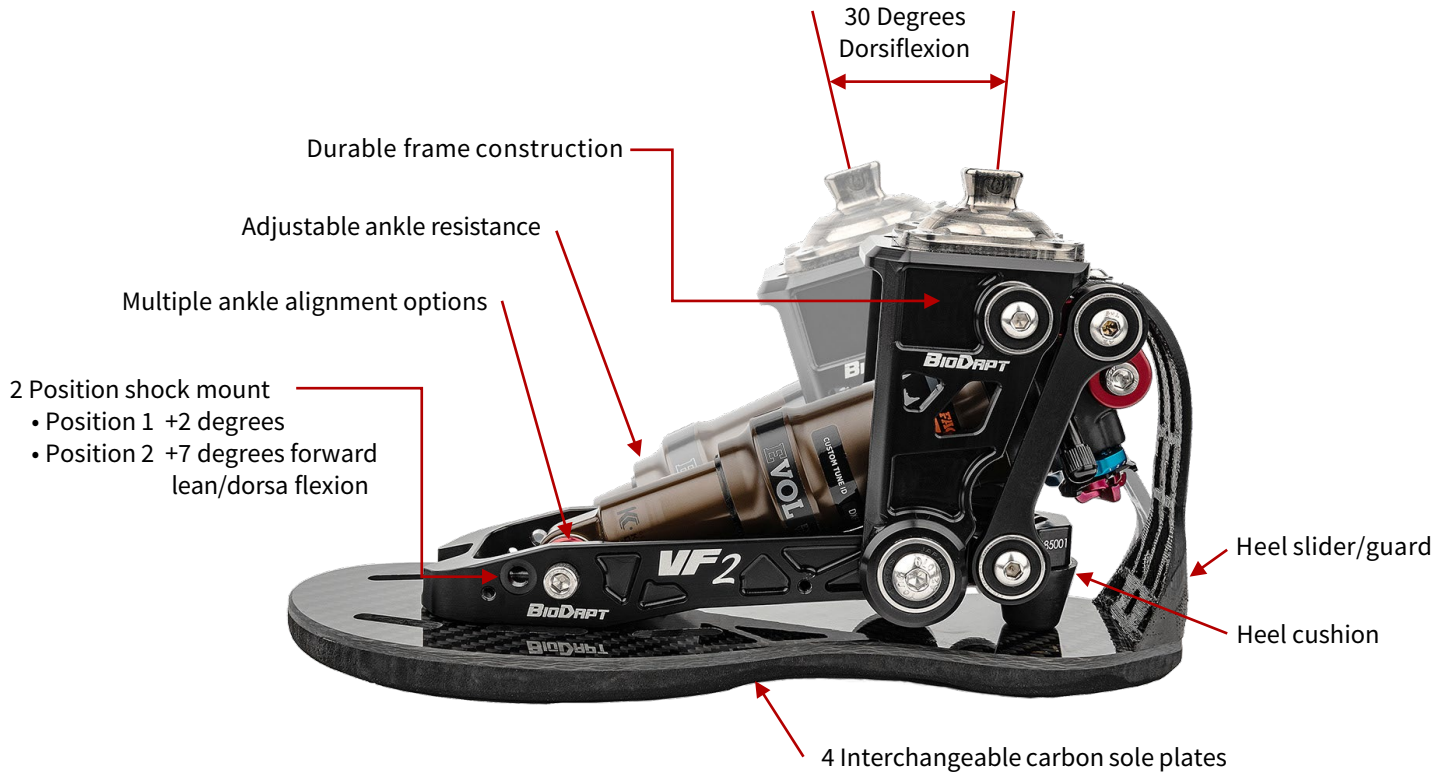

4 INTERCHANGEABLE SOLE PLATE OPTIONS

Get the most out of your sport or activity with versatile sole plate options.



FOAM SOLE

- General activity
- Used inside a shoe or boot



SHOE SOLE

- General activity
- Used without a shoe or boot



BIKE SOLE

- Cycling
- Used without a shoe or boot
- Cleat mount








MX SOLE

- Most durable boot sole
- Used without a shoe or boot

SPECIFICATIONS

SIZE RANGE	24 - 29 cm
EQUIPMENT WEIGHT	2 lbs 10 oz (1.19 kg)
MAX. PATIENT WEIGHT	225 lbs (102 kg)
BUILD HEIGHT	5.87 in - 6.125 in (14.92 - 15.56 cm)
FLEXION RANGE	30 degrees
KIT INCLUDES	Air pump, Multi-tool
WARRANTY	12 Months
WATERPROOF	Yes
WEATHERPROOF	Yes
ADAPTORS	Foot has proximal male pyramid
ACTIVITY LEVEL	Unlimited, K4+

VERSA FOOT2 PART NUMBERS WITH SOLE PLATE OPTIONS

PART NUMBER	DESCRIPTION	
VF85-FM**R (Right) VF85-FM**L (Left) <i>** Indicates size in cm</i>	Versa Foot2 with Foam Sole Plate - High density, smooth foam sole to be used inside a shoe or boot	 **Sizes: 25-29cm
VF85-SH**R (Right) VF85-SH**L (Left) <i>** Indicates size in cm</i>	Versa Foot2 with Shoe Sole Plate - Vibram shoe tread for use without a shoe. (Foot Guard recommended)	 **Sizes: 25,26,27, 29cm
VF85-MXR (Right) VF85-MXL (Left)	Versa Foot2 with MX Sole Plate - Performance foot with durable rubber, raised heel for use on motorcycle/atv foot peg. (Foot Guard recommended)	 ** One Sizes: 25cm
VF85-BKER (Right) VF85-BKEL (Left)	Versa Foot2 with Bike Sole Plate - Specific sole plate for cycling. Integrated cleat mount for use with “clipless” pedal system. (Foot Guard recommended) (Cleat not included)	 ** One Sizes: 25cm
VF85-Custom**R (Right) VF85-Custom**L (Left) <i>** Indicates size in cm</i>	Versa Foot2 with Custom Sole Plate - We can also build custom sole plates for specific applications such as Heavy Duty, raised heel, low profile etc.. Call us if you want something different than listed and we can look into building it.	
FG85000	Versa Foot2 Guard—protective cover.	
HLG85000	Versa Foot2 Heel Guard/Slider. Mounts to the heel end of sole plate. Included with all Sole Plate purchases	

1 Year Limited Warranty on foot frame and components

1 Year Limited Warranty on FOX SHOCK

INDIVIDUAL SOLE PLATE OPTIONS

PART NUMBER	DESCRIPTION	
SP85-FM**R (Right) SP85-FM**L (Left) <i>** Indicates size in cm</i>	Foam Sole Plate - High density, smooth foam sole to be used inside a shoe or boot - Heel Guard is included with sole plate (not shown)	 **Sizes: 25-29cm
SP85-SH**R (Right) SP85-SH**L (Left) <i>** Indicates size in cm</i>	Shoe Sole Plate - Vibram shoe tread for use without a shoe. - Heel Guard is included with sole plate (not shown) (Foot Guard recommended)	 **Sizes: 25,26,27, 29cm
SP85-MXR (Right) SP85-MXL (Left)	MX Sole Plate - Durable rubber, raised heel for use on motorcycle/atv foot peg. - Heel Guard is included with sole plate (not shown) (Foot Guard recommended)	 ** One Sizes: 25cm
SP85-BKER (Right) SP85-BKEL (Left)	Bike Sole Plate - Specific sole plate for cycling. Integrated cleat mount for use with “clipless” pedal system. Heel Guard is included with sole plate (not shown) (Foot Guard recommended) (Cleat not included)	 ** One Sizes: 25cm
SP85-Custom**R (Right) SP85-Custom**L (Left) <i>** Indicates size in cm</i>	Custom Sole Plate - We can also build custom sole plates for specific applications such as Heavy Duty, raised heel, low profile etc.. Call us if you want something different than listed and we can look into building it.	
501-4816	Versa Foot2 Replacement Heel Bumper SOFT - YELLOW	
501-4817	Versa Foot2 Replacement Heel Bumper MEDIUM - RED	
501-4818	Versa Foot2 Replacement Heel Bumper HARD - BLUE	



Contact us for BioDapt prosthetics information. We personally understand the challenges our customers face finding the right prosthetics and strive to make the process easy.



320.252.0894



BIODAPTINC.COM



5454 CTY 141
ST. CLOUD, MN 56301



INFO@BIODAPTINC.COM

“ I have two main passions in my life, being a competitor and designing/creating things in my shop. BioDapt allows me to do both while helping others perform at their highest level.

Mike Schultz, President BioDapt ”



OUR STORY STARTS HERE

The BioDapt story starts with Mike Schultz, a professional snowmobile racer that lost his leg above the knee to an injury while racing in 2008. Unable to find prosthetic equipment to get back in action and refusing to compromise, he created his own. His new designs allow him to perform at a high level in all the sports he loves including snowboarding, motocross, and snowmobile racing. BioDapt was created in 2010 to help other amputees achieve their goals.

MIKE SCHULTZ | President
PARALYMPIC GOLD AND SILVER MEDALIST (SNOWBOARDING)
10 TIME X GAMES GOLD MEDALIST (MOTOCROSS, SNOWMOBILE, SNOWBIKE)